

HARTFORD
BED RANGE



HF-136

4 Section Manual Bed

- Height Adjustable (Foot Operated)
- Backrest (Gas Spring Assisted)
- Knee Section (Gas Spring Assisted)
- Calf Section (By Ratchet Mechanism)
- Forward and Reverse Trendelenburg (Gas Spring Assisted)
- 125mm Central Locking Castors
- Integrated Safety Sides
- Removable ABS Head & Foot Boards.



HF-133

2 Section Manual Bed

- Fixed Height
- Backrest (Gas Spring Assisted)
- Forward and Reverse Trendelenburg (Gas Spring Assisted)
- 125mm Castors (2 x Individual Locking & 2 Free Swivel)
- Integrated Safety Sides
- Removable ABS Head & Foot Boards.



HF-135

2 Section Manual Bed

- Height Adjustable (Foot Operated)
- Backrest (Gas Spring Assisted)
- Forward and Reverse Trendelenburg (Gas Spring Assisted)
- 125mm Central Locking Castors
- Integrated Safety Sides
- Removable ABS Head & Foot Boards.



HF-132

4 Section Manual Bed

- Fixed Height
- Backrest (Gas Spring Assisted)
- Knee Section (Gas Spring Assisted)
- Calf Section (By Ratchet Mechanism)
- 125mm Castors (2 x Individual Locking & 2 Free Swivel)
- Integrated Safety Sides
- Removable ABS Head & Foot Boards.



HF-134

4 Section Manual Bed

- Fixed Height
- Backrest (Gas Spring Assisted)
- Knee Section (Gas Spring Assisted)
- Calf Section (By Ratchet Mechanism)
- Forward and Reverse Trendelenburg (Gas Spring Assisted)
- 125mm Castors (2 x Individual Locking & 2 Free Swivel)
- Integrated Safety Sides
- Removable ABS Head & Foot Boards.



HF-131

2 Section Manual Bed

- Fixed Height
- Backrest (Gas Spring Assisted)
- 125mm Castors (2 x Individual Locking & 2 Free Swivel)
- Integrated Safety Sides
- Removable ABS Head & Foot Boards.





Calf Section
This part of the 4 section bed can be raised and lowered manually via a ratchet mechanism, but only when the knee break section is raised. It automatically flattens when the knee section is levelled.



Safety Sides
¾ length integrated collapsible safety sides are fitted to all beds. In full compliance with the requirements of the international standard EN 60601-2-52.



Drainage Bag Hooks
On all models of the HARTFORD bed, a rail is provided on each side of the mattress platform complete with hooks to hold a drainage bag.




Head/Foot Boards
Lightweight ABS head and foot boards can be quickly and easily removed allowing access to the head and neck regions during general nursing and in the event of an emergency.




Trendelenburg
By gripping an ergonomically designed handle and pulling a release lever, the mattress platform can be tilted, with gas spring assistance, to any angle between +12° and -8°.



Accessory Fittings
Four sockets, two at each end of the bed, are provided to accommodate a lifting pole (head end only), IV poles and various traction and orthopaedic frames.

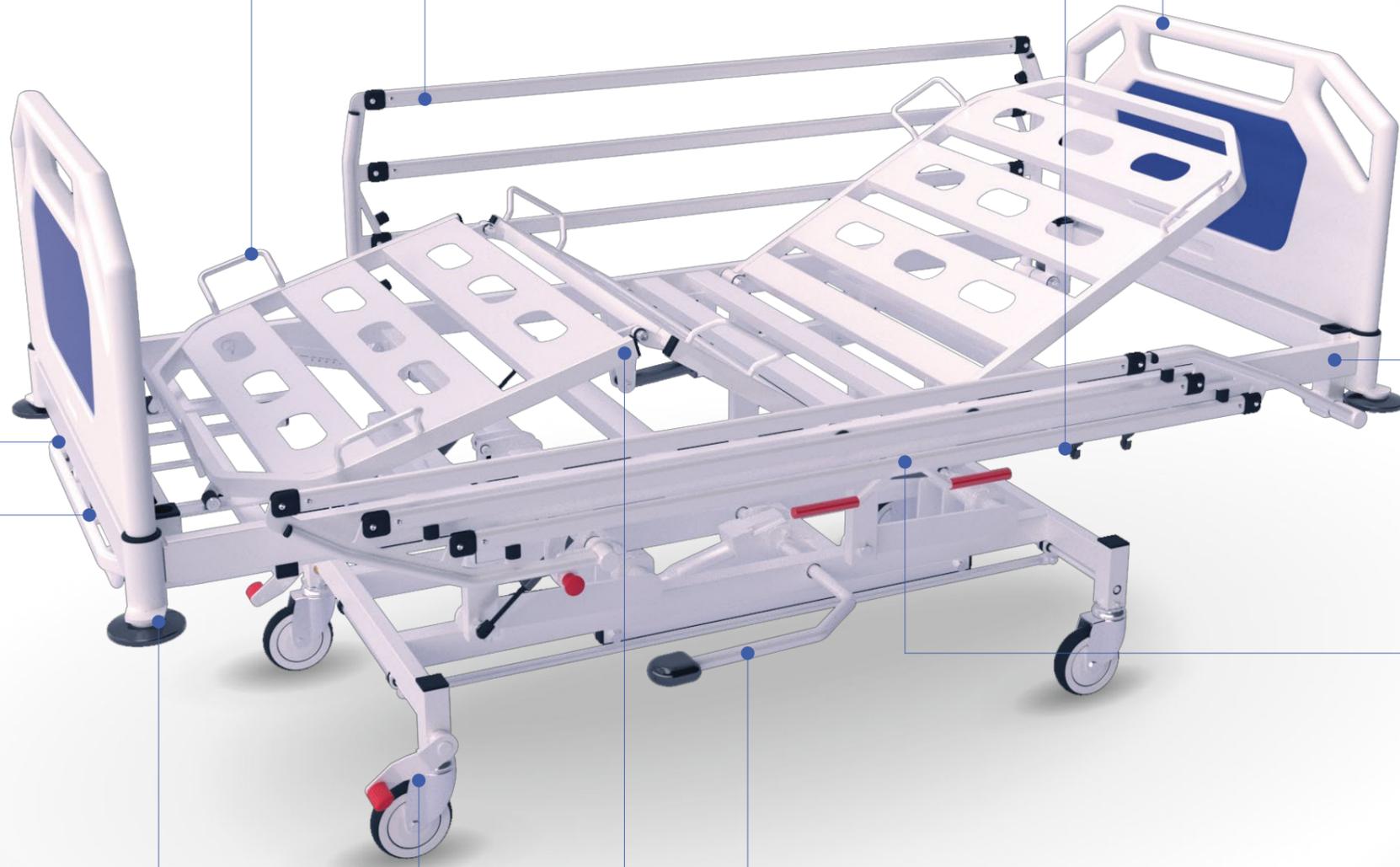

Linen Maker
A 'pull out' linen maker is provided on most HARTFORD models at the foot end of the bed to allow linen to be conveniently stored whilst the bed is being changed.

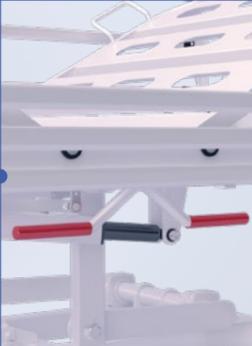


Backrest & Knee-Break
A handle on each side of the bed is provided to raise the gas spring assisted backrest and knee break.

Corner Bumpers
Non-marking 100mm diameter rotating corner bumpers are standard on all HARTFORD models. They protrude sufficiently to ensure that the bed is protected when being transported through narrow doorways.



Castors
Castors can be fully locked by downward activation of either brake pedal at the foot end of the bed. A steering lock can also be activated on one castor at the head end, making it easier to manoeuvre and steer the bed.

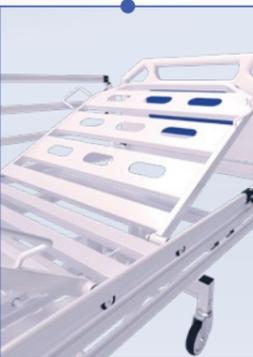

Juvenile Bed Option
The mattress platform for all Hartford bed models is constructed of sheet steel panels and are available in a juvenile size of 1800 x 760mm. All other bed features remain the same.



Height Adjustment
Adjustment of the height of the bed is achieved using a hydraulic foot pump. Pedals are located on each side of the bed which when repeatedly pressed effortlessly raise the bed. Lifting the pedal lowers the bed.

STANDARD FEATURES AND OPTIONS

ACCESSORIES

	BED MODEL					
	HF-136	HF-135	HF-134	HF-133	HF-132	HF-131
Juvenile Bed Size (1800 x 760mm)	○	○	○	○	○	○
125mm Central Locking Castors	✓	✓	○	○	—	—
150mm Central Locking Castors	○	○	○	○	—	—
125mm Swivel & Brake Castors	○	○	✓	✓	✓	✓
Pull Out Linen Maker	✓	✓	✓	✓	○	○
Bed Extension	○	○	○	○	○	○
Sheet Steel Slats for Mattress	✓	✓	✓	✓	✓	✓
Wire Mesh Mattress Platform	○	○	○	○	○	○
Urine Bag Holder	✓	✓	✓	✓	✓	✓
Colour - Grey White RAL 9002	✓	✓	✓	✓	✓	✓

✓ Standard
 ○ Optional
 — Not Available

TECHNICAL SPECIFICATIONS

HF-131	90kg	Overall Dimensions	2145 x 1035mm	
HF-132	98kg	Mattress Platform (Adult)	2020 x 900mm	
HF-133	108kg	Mattress Platform (Juvenile)	1800 x 760mm	
HF-134	116kg	Variable Height Models	460 to 760mm	
HF-135	125kg	Fixed Height Models	680mm	
HF-136	133kg	Adult Mattress (L x W)	2000 x 900mm	
		Juvenile Mattress (L x W)	1790 x 760mm	
		Mattress Thickness (Max/Min)	150/100mm	
Backrest Angle	0 - 70°		Safe Working Load	220kg
Knee Break Angle	0 - 35°		Maximum Patient Weight	180kg
Trendelenburg Angle	0 - 12°			
Reverse Trendelenburg	0 - 8°			

Hospi-Furn reserves the right to change technical specifications without prior notice.



HA-040 | Full Orthopaedic Frame

HA-050 | Small Traction Frame




HA-021 | 2 Hook IV Pole




HA-020 | 4 Hook IV Pole



HA-010 | Lifting Pole

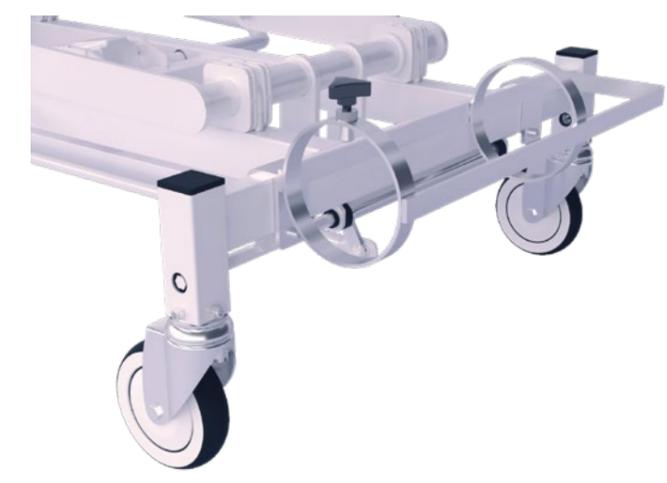

HA-100 | Oxygen Bottle Holder

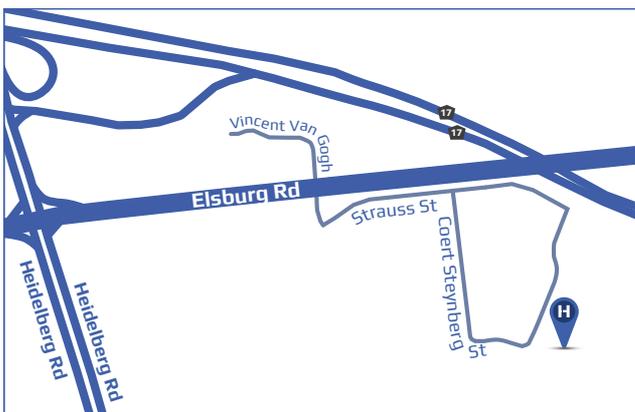

-26.261264, 28.332704

Physical Address:

23 Coert Steynberg St, Van Eck Business Park, Brakpan, Gauteng

Postal Address:

P.O. Box 6920, Westwood, Gauteng, 1477

Contact Details:

Tel: 011 894 1600

Fax: 086 557 1577

Email: sales@hospi-furn.co.za

Web: www.hospi-furn.co.za